

# Infrared Forehead Thermometer

4DET-306

## Infrared forehead thermometer designed for measuring body temperature

The 4DET-306 infrared forehead thermometer is specifically designed for use on the forehead. It is a no-contact device capable of measuring people's body temperature by detecting the intensity of infrared light emitted from the forehead. It converts the measured heat into a temperature reading displayed on the LCD.

### Technology Highlights

- Fast reading (approx. 1 second)
- No probe covers needed
- Beeps when reading is complete
- Mercury free
- Stores 10 readings
- Automatic Shut-off


**No-Contact Device!**



### 4DET-306 Manufacturer's Specifications

<b>Temperature Range:</b> Infrared: 89.6° to 109.4°F (32° to 43°C) Object: 32° to 212°F (0° to 100°C)
<b>Accuracy:</b> Forehead Mode: $\pm 0.4^{\circ}\text{F}$ ( $\pm 0.2^{\circ}\text{C}$ ) between 95.9°F to 107.6°F (35.5°C to 42°C) at 59°F to 95°F (15°C to 35°C). Object Mode: $\pm 4\%$ or $\pm 0.4^{\circ}\text{F}$ ( $\pm 0.2^{\circ}\text{C}$ ) whichever is greater.
<b>Resolution:</b> 0.1°F/0.1°C
<b>Response Time:</b> 1 second (approximately)
<b>Power:</b> (2) AAA DC3V battery
<b>Battery Life:</b> 6000 readings (approximately 1 year)
<b>Dimensions (L x W x D):</b> 6" x 1.6" x 1.7" (153 mm X 41 mm X 44 mm)
<b>Weight (with batteries):</b> 3 oz (84 grams)
<b>Warranty:</b> 1 year
<b>UPC Code:</b> 6970392215631

FDA  
510K Cleared

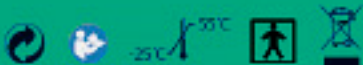

**Manufactured by JoyTech Healthcare Co., Ltd.**  
No.365, Wuzhou Road,  
Yuhang Economic Development Zone,  
311100, Hangzhou, China  
T +86-571-81957767

**Distributed by Cooper-Atkins**  
Emerson Commercial & Residential Solutions  
33 Reeds Gap Road | Middlefield, CT 06455  
T +800-835-5011 | F +860-347-5135  
www.cooper-atkins.com | 67-1791 | 0420

Made in China



Infrared Forehead Thermometer measures core body temperature and is designed to measure the temperature of the skin surface over the temporal artery, a major artery of the head. Therefore, body temperature changes are generally reflected sooner in the forehead than they are in other parts of the body such as oral, rectal and underarm.

For additional information contact your local foodservice distributor



CONSIDER IT SOLVED.™